

SmartRelease®

Endoscopic Soft Tissue Release System

For Carpal Tunnel Release



MICROAIRE®
For Surgery. For Life.™

Minimally Invasive Endoscopic Treatment for Carpal Tunnel Syndrome

SmartRelease[®]

Endoscopic Soft Tissue Release System

- Easy to learn, simple to use
- Documented safety and efficacy ^{1,2,3}
- Complete set of custom instrumentation
- Surgical protocol developed and tested by a panel of hand surgeons
- Designed exclusively for transecting the transverse carpal ligament



Standard Disposable Blade Assembly

(Single) | 81010-1
(6-pack) | 81010-6

- Sterile package
- Window at the tip of the assembly provides direct view of the transverse carpal ligament and surrounding anatomy
- Disposable blade assembly features a retractable blade that rises straight up, without curving outward



Onyx™ Low Profile Disposable Blade Assembly

(Single) | 83030-1
(6-pack) | 83030-6

- Less volume for smaller, tighter carpal tunnels
- Sharper blade than the standard blade assembly
- Scraper tip assists with clearing the tenosynovium
- Black color absorbs light for improved tissue contrast





SmartRelease[®] Aluminum Handpiece | 83014

- Aluminum, pistol-grip handpiece
- Endoscope slides into position from the back
- Blade assembly click-locks onto the front
- Turn-latch on top of the instrument allows blade assembly to be rotated for right- or left-hand surgery and locks the blade in any position along the rotation axis



2.9mm Eyepiece Autoclavable Endoscope | 81025

- Connects to a standard camera coupler
- Storz[®] and Wolf[®] light-post adapters

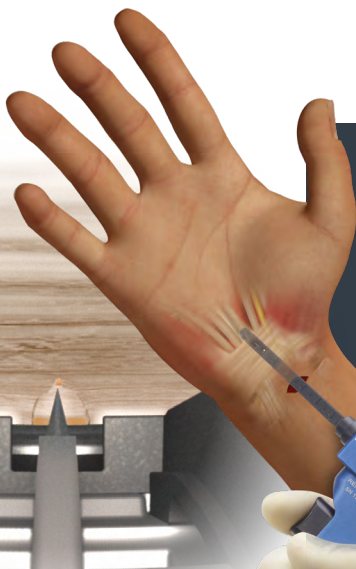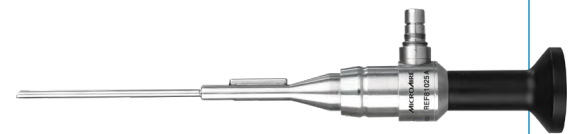

“ Single-portal endoscopic surgery is a safe and effective method of treating carpal tunnel syndrome.²