

Coaxial

Micro

Stick[®]

Introducer Set

Performance
By Design

Coaxial Micro-Stick®


Medcomp® introduces our new

Through advanced engineering and careful design details, we focused on the system as a whole to create performance by design.

Micro-Stick® Design

Engineers designed the tip around enhancing the performance of the procedure.

- Longer tip for a more gradual dilation of the vessel.
- New, sleeker tip design features a zero-clearance transition between sheath and dilator, allowing a smoother insertion into the vessel.
- Stiff version features a steel cannula inside the dilator, providing extra rigidity for accessing tough vessels.

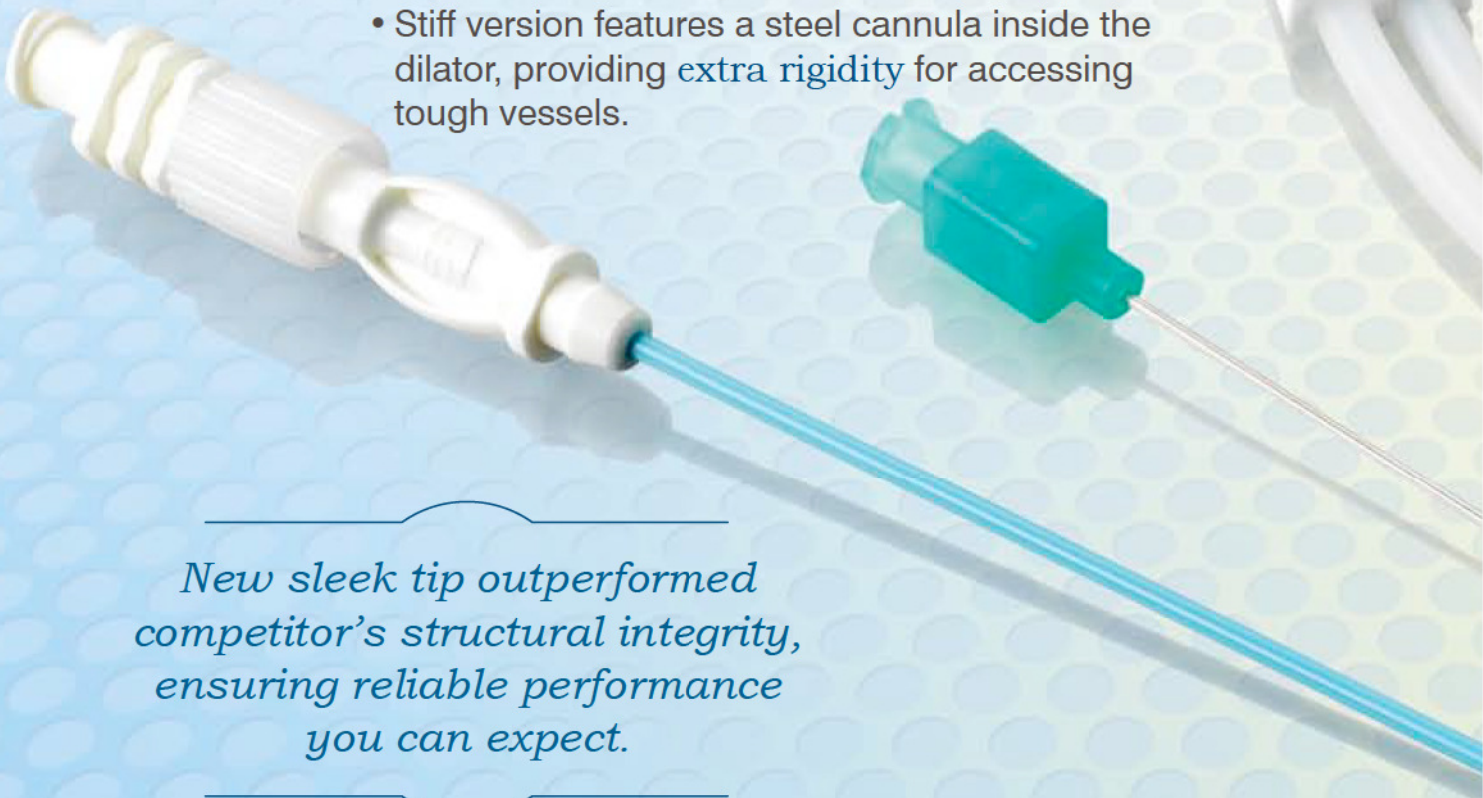

New sleek tip outperformed competitor's structural integrity, ensuring reliable performance you can expect.

New
STIFF Version
Now Available.

new coaxial Micro-Stick® design.

Medcomp® offers several options, so clinicians can easily tailor their procedure to the patients' needs.

Coaxial Micro-Stick®

- A choice of 4 and 5 French sheath sizes.
- Stiffer version available to ensure access into tough vessels.

21 Gauge Super Sharp™ Needle

- Provides less traumatic insertion.
- Echogenic tip for higher visibility.
- Bevel indicator for confidence in bevel position.
- Echogenic and non-echogenic tips available.

Variety of .018 x 45cm Guide Wires

- Nitinol or Steel wires.
- Radiopaque distal tip enhances visualization.
- Curvature of floppy tip allows flexible ease into the vessel.